


Don't just seal, SURESEAL

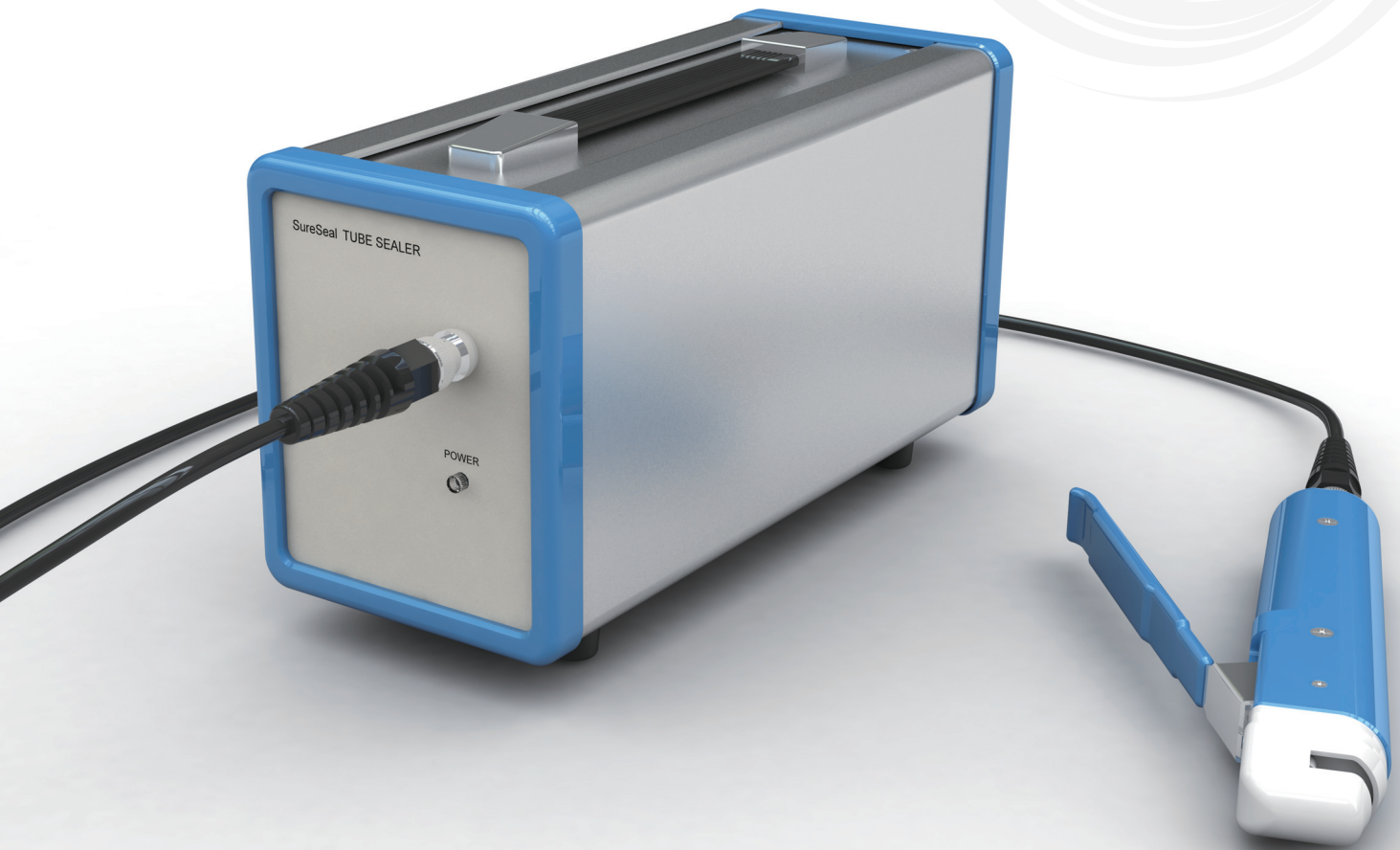

CENTRON

***SureSeal*[™] SE450**

TUBE SEALER

All-round hand unit tube sealer

SureSeal[™]SE450 is a compact and lightweight tube sealer with a hand unit. The machine is suitable for use where work space is limited or where one tube sealer machine will be shared by multiple users

■ Features

- AC-powered, with handheld sealing unit
- Ergonomic design with standard 1.8 meters (6 ft) cable
- Seals wide range of blood bag tubing - good for tubing diameter 2-6 mm (1/4")
- Fast sealing in one second / No warm-up required
- Makes seals with a notch at the center for easy segmenting of tubing
- Head protector is detachable for sealing head cleaning
- Suitable at blood collection sites or laboratories where multiple users should share a machine

Centron

Smart Solutions for Smarter Bloodbanking

Main unit

- The machine is suitable for use where work space is limited or where one tube sealer machine will be shared by multiple users. When the trigger lever is pulled, the machine generates strong high-frequency power to melt thermoplastic tubing. It makes a seal on tubing in 1-2 seconds.



RF cable connector


- Connect the RF cable between this connector and the same connector on the hand unit.

Power indication lamp

- The green LED lamp lights when the power switch is on.

Sealing hand unit

- Additional plug-in type *SureSeal*™ sealing hand unit
- Ergonomic design with standard 1.8m (6ft) cable
- Satisfying various operating conditions

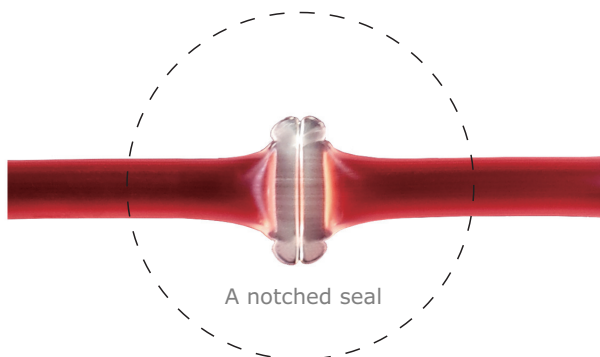

Specification

Product	TUBE SEALER
Model	SE450
Power	220-240V~ 50/60Hz 250W
Protection	FUSE, 2A/250V~
Dimension	Main unit: 106(W)mm/206(L)mm/150(H)mm
	Head unit: 28(W)mm/204(L)mm/44(H)mm

Information



DISTRIBUTOR

Centron Centron Corporation

Address: 44 Beolmal-ro, Dongan-gu, Anyang, Gyeonggi-do,
Postal Code 431-804 / Republic of Korea
Phone: +82-31-422-7890 / **Fax:** +82-31-422-7866
E-mail: info@centrons.com
Website: www.centrons.com



※ The specification and design are subject to change without prior notice for product upgrades.